iQ-SR2PDF CONVERTER 1.2

From SR to PDF

Automate the creation of your radiology reports

DICOM Structured Reports have significantly improved medical diagnostic report communication as they are readable by medical information systems worldwide and lead to a much better report turnaround time. But outside of the radiology department, referring and treating physicians still rather prefer to receive the results for their patient's procedures in a PDF or paper format.

iQ-SR2PDF CONVERTER helps to combine the efficiency of a DICOM SR with the traditional expectations of receiving a radiology report in letter form.

This smart application complements our iQ-SYSTEM PACS product range and offers iQ-4VIEW and iQ-VIEW users further advantages as it automatically transforms DICOM Structured Reports from both applications into customized PDF letters and makes those available for sharing.

BENEFITS



Save valuable time and automatically transform iQ-4VIEW and iQ-VIEW DICOM SR into DICOM PDF



Share the reports with multiple DICOM targets (using iQ-WEB)



Create your report templates per imaging center (institution), radiology reading group (reading entity) or radiologist (reading physician)



Digitally sign the reports by adding single (verifier) or dual signature (verifier and observer)




Manage the PDF reports for referring physicians in iQ-WEB or any other PACS supporting DICOM PDF objects



Use patient ID barcodes (EAN-13 and UPC-A supported)

How does iQ-SR2PDF CONVERTER fit into the medical imaging workflow?

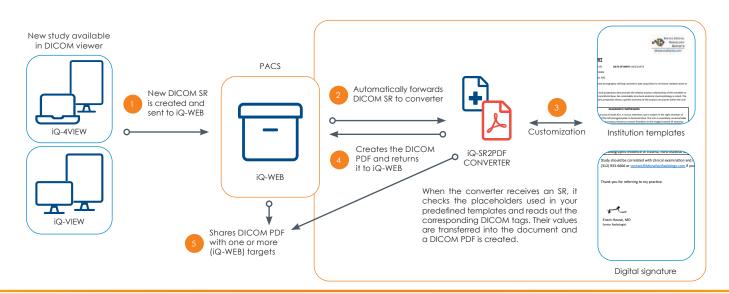

iQ-SR2PDF CONVERTER 1.2

SUPPORTED PLACEHOLDERS

Placeholders can be used to automatically include information provided in the DICOM SR in the radiology report. Every placeholder matches a specific DICOM tag. The value in this DICOM tag is read out and added to the letter text in the position where the placeholder is set.

The following placeholders are available*:

- Patient demographics (patient name, gender, date of birth and patient ID), patient's address and phone number
- Referring physician's name, address and phone number, requesting physician's name, admitting diagnoses description
- Study description, study date and time, accession number, body parts examined, modalities in study, series description list, institution name, device manufacturers and model names
- Report status and verification date, report text fields (Technique, Observation, Summary), observer's and verifier's names, incl. signatures and professional roles
- Barcode of patient ID (EAN-13 or UPC-A)
- * iQ-4VIEW ≥ 2.2.3 is required for full use of all placeholder options. The use of placeholders is limited for iQ-VIEW.

PDF LETTER EXAMPLES







General Requirements: Installation on any currently supported 64 bit Microsoft Windows OS; Microsoft Word ≥ 2007 (for letter templates and appropriate conversion); DICOM SRs created by iQ-4VIEW and iQ-VIEW (unsupported reports are rejected)

Stamp of sales partner